

# SprintRay Dental Model

## Instructions For Use

### Indications for Use

The SprintRay Dental Model is a light-curable polymerizable resin intended to be used for fabricating teeth models, which can be used for multiple purposes such as:

- Study models
- Removable die model
- Models to fabricate aligners
- Models to fabricate retainers
- Models to fabricate bleaching trays
- Wax-up models

This material is an alternative to traditional dental model material.

---

### Contraindications

SprintRay Dental Model is contraindicated when:

- a patient is known to be allergic to any of the ingredients
- there is direct intraoral contact with resin that is not fully cured
- it is used for any purpose other than its indications for use

---

### Device Description

Dental Model is an alternative to traditional material for the fabrication of dental models. It is intended exclusively for professional dental work.


---

### Printing and Hardware Parameters

These device specifications have been validated using the following manufacturing products. Any products or processes not specified in this document are outside of the device specifications.

- CAD File:** CAD file of treatment device in STL file format
- Printer:** SprintRay Pro or Pro S or Pro 2 3D printer
  - Pro or Pro S: 55 or 95-micron XY resolution
  - Pro 2: 35 micron XY resolution
- Software:** RayWare Desktop or RayWare Cloud
  - STL file import
  - Manual/automatic orientation
- Printing Parameters**

- i. For optimized accuracy, orient the intaglio surface facing away from build platform as shown in the photo below
- ii. Select the desired layer thickness (RayWare will typically default to 100 microns)
- iii. Default support structures

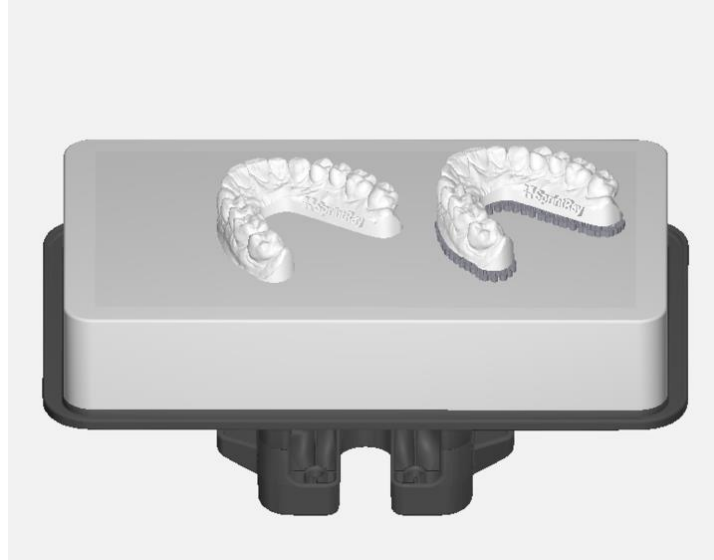

- e. **Wash Device:** SprintRay ProWash S or SprintRay Pro Wash/Dry
  - i. 91% or higher IPA
  - ii. Standard preprogrammed multi-cycle wash
- f. **Cure Device:** SprintRay ProCure 2 or NanoCure
  - i. Use SprintRay-recommended curing times that are built in the device

---

## Warning and Precautions

SprintRay Dental Model is non-toxic in processed, cured form, and is classified as a biocompatible material. In uncured form, Dental Model is classified as a sensitizer. When washing with solvent or grinding the device, do so in a well-ventilated area with proper protective equipment. Wear protective gloves, clothing, eyewear, and face protection when handling.

- **Skin Contact:** May cause skin irritation. If unprocessed resin contacts skin, wash thoroughly with soap and water. May cause an allergic skin reaction. If skin sensitization occurs, stop using. If dermatitis or other symptoms persist, seek medical assistance.
- **Inhalation:** High vapor concentration may cause headache, irritation of eyes and/or respiratory system. If exposed to a high concentration of vapor or mist, move to fresh air. Use oxygen or artificial respiration as required.
- **Eye Contact:** Wash the contacted area thoroughly with soap and water.

- **Ingestion:** Contact your regional poison control center immediately.
- **Use of Incompatible Components:** Do not substitute any of the components of the device system, i.e., device photopolymer materials, bonding systems, scanners, 3D printers, post-curing units, CAD/CAM software, templates, and tools. Use only those specifically identified in this labeling. Unauthorized changes may result in a device that is outside of specification. Contact the manufacturer for compatible components.
- Maintain and calibrate equipment according to manufacturer instructions.
- **Minor Color Differences:** Shade variance may occur due to inadequate shaking and mixing of the original packaging before use; inadequate stirring in the resin tank before use; insufficient post-curing

---

## Storage

- **Material Reuse:** The remaining resin in the resin tank can be reused. You may use a filter to ensure the resin is free from any cured particles to avoid print failures. The remaining material in the tank can be poured back into the resin bottle upon filtration. This process can be repeated until the material in the bottle is fully consumed. Please note that in the case of reuse, the resin must be filtered and poured back into the same bottle.
- Store Dental Model at 15-25°C (60-77°F) and avoid direct sunlight.
- Keep the bottle closed and/or the tank lid securely attached when not in use.
- Do not use Dental Model after the expiration date printed on the bottle.
- Resin must be protected from exposure to light, as spontaneous polymerization is possible. The bottle must be tightly closed after every usage.



Do not use expired resin as biocompatibility, performance, and print stability may be compromised.

---

## Fabrication of Device

This resin was validated using the following workflow. Failure to follow the recommended practices may lead to undesired safety and performance implications.

Any deviation from these instructions for use may negatively affect the physical and/or chemical qualities of the resin and the biocompatibility of the final device.

If applicable, refer to the Workflow Guide for detailed best practices for producing specific appliance types with SprintRay resins.

## Designing

The device is designed in STL file format by a dental design service, preferably SprintRay Cloud Design, or other dental CAD software using digital anatomical data from the patient. This STL file is delivered to the clinician for fabrication.

### **3D Printing**

Sign in to RayWare Cloud and select the appliance type; the algorithm will automatically orient and add supports. Select this material and use the desired layer thickness. Queue the job to your printer.

Shake the resin bottle thoroughly for one minute, then pour into the resin tank up to at least the min fill line. From the printer touchscreen, assign the resin tank to the proper material and shade, then navigate to the printer queue. Start the print job.

### **Part and Support Removal**

After your device has been printed, remove it from the print platform using the provided Print Removal Tool. Remove all supports using a flush cutter or round diamond disc. Cut as close as possible to the device to minimize the smoothing and finishing procedure.

### **Washing and Drying**

- Use  $\geq 91\%$  IPA to wash the device using the SprintRay ProWash S or SprintRay Pro Wash/Dry: Standard cleaning multi-cycle

To ensure the proper function of the wash unit, always follow on-screen instructions for device cleanliness and maintenance. Dry the part completely before post-curing.

### **Post Curing**

Use one of the following post-curing equipment from SprintRay to cure the device and select the preprogrammed profile for Dental Model:

- ProCure 2 (preprogrammed material profile)
- NanoCure (preprogrammed material profile)

### **Finishing**

Wash and clean the device in an IPA sonication bath for 2 minutes following with a brush using soap and warm water.

---






















## **Disposal Considerations**

Always follow federal, state, and local regulations for hazardous waste disposal. To ensure proper classification, consult your local regulations. US guidelines can be found in 40 CFR part 261.3. Liquid resin must be cured completely before regular disposal. Simply pour it into a clear

container and expose it to direct sunlight until hardened or in one of the post-cure boxes. SprintRay Dental Model is not an environmental hazard in its final, fully cured state. Once cured, it can be thrown away with regular trash.

## Symbol Guide

The below table provides reference for symbols that may appear on the resin bottle label.

	Keep away from sunlight		Use-by date
	Consult instructions for use		European conformity
	LOT number		SKU number
	Manufacturer		Temperature limit
	Prescription only		Medical device
	Environmental hazard		Irritation
	Unique device identifier		Importer
	Indicates the authorized representative in Switzerland		Authorized representative in the European community
	Manufacturing date		Wear gloves
	Health hazard		UK Conformity Assessed (UKCA) Marking
	UK responsible person		

## Additional Help & Support

We are here to support you throughout the implementation period of your new technology. Our experienced support technicians are available M - F from 6 AM - 5 PM PT at 800-914-8004.

---

## Contact Information

For product assistance, please review help information at: <https://sprintray.com/digital-dentistry/>

To report product issues, please contact SprintRay at: <https://support.sprintray.com/hc/en-us/requests/new>

Phone: 1-800-914-8004

---



## Manufacturer information

SprintRay Inc.  
2710 Media Center Dr., Suite #100A  
Los Angeles, CA 90065, USA